

Naída S

UltraPower, SuperPower and CRT hearing aids

User Guide



PHONAK
life is on

Contents

1. Welcome	5
2. Important safety information: Please read before operating your hearing aid for the first time	6
Hazard warnings	6
Information on product safety	9
3. Hearing aid description	12
4. Step-by-step instructions for using the hearing aid	16
Step 1. Getting to know your hearing aid	16
Step 2. Battery preparation	19
Step 3. Inserting the battery	20
Step 4. Turning on the hearing aid	21
Step 5. Inserting the hearing aid	22
Step 6. Adjusting the volume control behind the ear	28
Step 7. Selecting the hearing program behind the ear	30
Step 8. Removing the hearing aid	32
Step 9. Turning off the hearing aid	34

Step 10. Changing the battery	35
Step 11. Changing the battery in hearing aid with a secured battery compartment	36
5. Care and maintenance	42
6. Troubleshooting	46
7. Wireless accessories and FM systems	48
Wireless accessories	48
FM systems	49
FM system description	51
Design-integrated FM receiver	51
Audio shoe with universal FM receiver	53
8. Service and warranty	57
9. Compliance information	59
10. Information and explanation of symbols	62

This user guide is valid for:

Wireless models	CE mark applied
Naída S IX SP	2011
Naída S IX UP	2011
Naída S V SP	2011
Naída S V UP	2011
Naída S III SP	2011
Naída S III UP	2011
Naída S CRT IX SP	2011
Naída S CRT IX UP	2011
Naída S CRT V SP	2011
Naída S CRT V UP	2011
Naída S CRT III SP	2011
Naída S CRT III UP	2011

Non wireless models	CE mark applied
Naída S I SP	2012
Naída S I UP	2012

1. Welcome

Your new hearing aid is a premium Swiss quality product. It was developed by Phonak, one of the world's leading companies in hearing technology. Your hearing aid offers the most modern digital hearing technology currently available on the market.

Please read these user instructions thoroughly in order to benefit from all the possibilities your new hearing aid has to offer. With proper use and care, you will be able to enjoy your hearing aid for many years to come.

If you have any questions, please consult your hearing care professional.

Phonak – life is on

www.phonak.com




2. Important safety information

Please read the information on the following pages before using your hearing aid.

A hearing aid will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions. Infrequent use of a hearing aid does not permit a user to attain full benefit from it. The use of a hearing aid is only part of hearing habilitation and may need to be supplemented by auditory training and instruction in lipreading.

2.1 Hazard warnings

 The intended use of hearing aids is to amplify and transmit sound to the ears and hereby compensate for impaired hearing. The hearing aids (specially programmed for each hearing loss) must only be used by the intended person. They should not be used by any other person as they could damage hearing.

-
- ⚠ Changes or modifications to the hearing aid that were not explicitly approved by Phonak are not permitted. Such changes may damage your ear or the hearing aid.

 - ⚠ Do not use the hearing aids in explosive areas (mines or industrial areas with danger of explosions, oxygen rich environments or areas where flammable anesthetics are handled).

 - ⚠ Hearing aid batteries are toxic if they are swallowed! Keep out of reach of children and mentally challenged people or pets. If batteries are swallowed, consult your physician immediately!

 - ⚠ If you feel pain in or behind your ear, if it is inflamed or if skin irritation and accelerated accumulations of ear wax occur, please check with your hearing care professional or physician.

 - ⚠ In very rare cases, the dome can remain in your ear canal when removing the hearing tube from the ear. In the unlikely case that the dome does get stuck in your ear canal, it is strongly recommended to see a physician for safe removal.

2. Important safety information

- ⚠ Hearing programs in the directional microphone mode reduce background noises. Please be aware that warning signals or noises coming from behind, e.g. cars, are partially or entirely suppressed.

- ⚠ The Hearing aid contains small parts that can cause choking, if swallowed by children. Keep out of reach of children and mentally challenged people or pets. If swallowed, consult a physician or hospital immediately.

- ⚠ External devices may only be connected if they have been tested in accordance with corresponding IECXXXX standards. Only use accessories approved by Phonak AG, to avoid any electrical shock.

- ⚠ The following is only applicable for persons with active implantable medical devices (i.e. pacemakers, defibrillators, etc.):
 - Keep the wireless hearing aid at least 15 cm (6 inches) away from the active implant. If you experience any interference, do not use the wireless hearing aids and contact the manufacturer of the active implant.

Please, note that interference can also be caused by power lines, electrostatic discharge, airport metal detectors etc.

- Keep magnets (i.e. battery handling tool, EasyPhone magnet, etc.) at least 15 cm (6 inches) away from the active implant.

2.2 Information on product safety

- ① Your hearing aid is resistant to water, sweat and dust under the conditions listed in chapter 5. After exposure to any of these elements the hearing aid must be cleaned and dried, as the hearing aid contains sensitive electronic parts.
- ① Protect your hearing aid from heat (never leave it near a window or in the car). Never use a microwave or other heating devices to dry your hearing aid. Ask your hearing care professional about suitable drying methods.

2. Important safety information

- ① When you are not using your hearing aid, leave the battery compartment open so that any moisture can evaporate. Make sure that you always completely dry your hearing aid after use. Store the hearing aid in a safe, dry and clean place.
- ① Do not drop your hearing aid! Dropping onto a hard surface can damage your hearing aid.
- ① Always use new batteries for your hearing aid. In case a battery is leaking, replace it immediately with a new one to avoid any skin irritation. You can return used batteries to your hearing care professional.
- ① Remove the battery if you are not using your hearing aid for a long period of time.

① Special medical or dental examination including radiation described below, may adversely affect the correct functioning of your hearing aids. Remove and keep them outside the examination room/area before undergoing:

- Medical or dental examination with X-ray (also CT scan)
- Medical examinations with MRI/NMRI scans, generating magnetic fields

Hearing aids don't need to be removed when passing security gates (airports etc.). If X-ray is used at all, it will be in very low doses, and will not affect the hearing aids.

① Do not use your hearing aid in areas where electronic equipment is prohibited.

3. Hearing aid description

The following user instructions describe three hearing aid models: The **UltraPower (UP)**, **SuperPower (SP)** and the **Canal Receiver Technology (CRT)** hearing aid.

The different hearing aids can be adapted with a variety of earmolds.


On the following pages use the labeled drawings to identify the hearing aid and the earmold that you use.

The following table shows what battery you need for your hearing aid model.

Model	Zinc air battery size (marking on package)	IEC code	ANSI code
SP	13 (orange)	PR48	7000ZD
UP	675 (blue)	PR44	7003ZD
CRT	13 (orange)	PR48	7000ZD

UltraPower and SuperPower model

Option A: with classic earmold

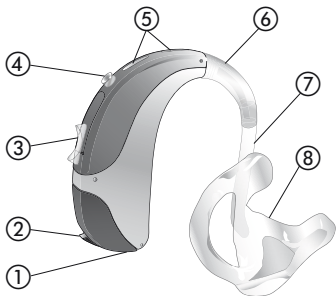


- ① Right or left indicator (right=red, left=blue)
- ② Battery compartment with ON/OFF switch
- ③ Volume control
- ④ Program button
- ⑤ Microphone input with Wind and Weather Protector
- ⑥ Hook/sound output
- ⑦ Hearing tube
- ⑧ Classic earmold (detachable)

3. Hearing aid description

UltraPower and SuperPower model


Option B: Power SlimTube with classic earmold

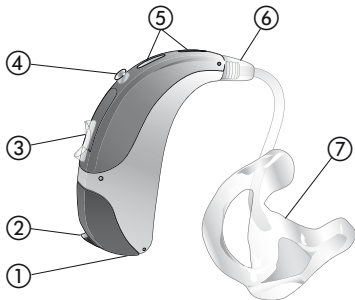


- ① Right or left indicator (right=red, left=blue)
- ② Battery compartment with ON / OFF switch
- ③ Volume control
- ④ Program button
- ⑤ Microphone input with Wind and Weather Protector
- ⑥ Power SlimTube
- ⑦ Classic earmold (detachable)

CRT model

Option A: with cShell

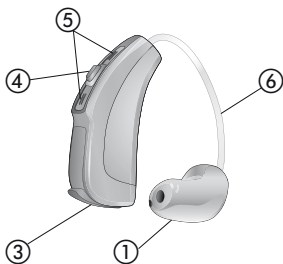

Option B: with Dome

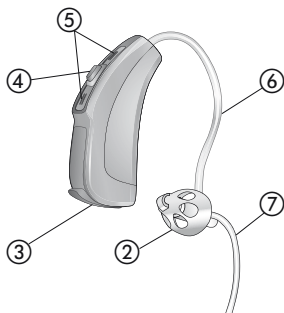

- ① cShell
- ② Dome (detachable)
- ③ Battery compartment with ON/OFF switch
- ④ Push button
- ⑤ Microphone input with microphone protection
- ⑥ Sound outlet
- ⑦ Retention

4. Step-by-step instructions for using the hearing aid

The following chapter provides step-by-step instructions on how to use your hearing aid. Please follow these steps carefully. In the drawings, the parts that are particularly important for each step are shown in green.

Step 1. Getting to know your hearing aid

Take the hearing aid in your hand and try out the hearing aid controls. This will make it easier for you to feel and operate the controls later on when wearing the hearing aid.

UltraPower and SuperPower model

To increase the volume, press the volume control upwards. To decrease the volume, press the volume control downwards (*Fig. 1a*).

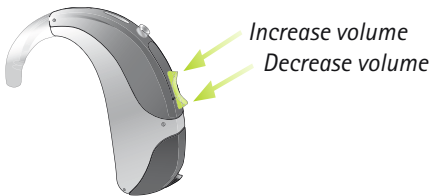


Fig. 1a

In order to change a hearing program, press the program button displayed (*Fig. 1b*).



Fig. 1b

4. Using the hearing aid

CRT model

Your hearing care professional can configure the push button as either a volume control or a program selector.

To increase the volume, press the button on your right hearing aid. To decrease the volume, press the button on your left hearing aid (*Fig. 1c*).

*Decrease
volume*



Left hearing aid

*Increase
volume*



Right hearing aid

Fig. 1c

In order to change a hearing program, press the button displayed (Fig. 1d).

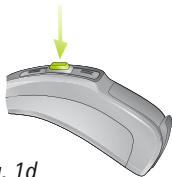


Fig. 1d

Step 2. Battery preparation

Remove the protective foil from the new battery (Fig. 2). Wait for two minutes before inserting to activate the battery.

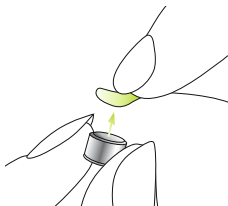


Fig. 2

4. Using the hearing aid

Step 3. Inserting the battery

Take the hearing aid in your hand (*Fig. 3a/3c*) and open the battery door. Insert the battery so that you see the + symbol on the battery (*Fig. 3b/3d*).

UP and SP model



Fig. 3a



Fig. 3b

CRT model



Fig. 3c



Fig. 3d

Step 4. Turning on the hearing aid

Switch on the hearing aid by closing the battery compartment (*Fig. 4a/b*).

UP and SP model

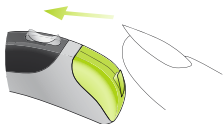


Fig. 4a

CRT model

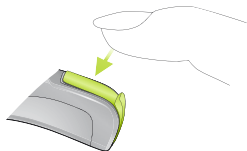


Fig. 4b

- ① Open and close the battery compartment gently and carefully.
- ① Check if you feel resistance when closing the battery compartment. Check to make sure the battery is correctly inserted and in the correct direction. If the battery is inserted the wrong way, the hearing aid will not work and the battery compartment can be damaged.

4. Using the hearing aid

The hearing aid is now switched on. It can take up to 15 seconds before it starts.

(When you wear the hearing aid, you will hear an acoustic signal).

Step 5. Inserting the hearing aid

This step shows you how to correctly insert the hearing aid. The instructions apply for the **UltraPower**, **SuperPower** and **CRT** hearing aids. The type of earmold determines how the hearing aid is inserted.

- ① Each hearing aid is individually programmed for your right or left ear. For the first step identify right and left so you will be able to insert your hearing aid in the correct ear.

Step 5.1 Identifying the hearing aid for the right and left ear

Color indicator:

right = red

left = blue

UP and SP model

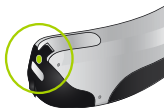


Fig. 5a

CRT model

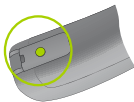


Fig. 5b

4. Using the hearing aid

Step 5.2 Inserting an UltraPower or SuperPower model with classic earmold

Before inserting the earmold, check that the earmold tube is slid over the hearing aid hook correctly. With your right hand, take your earmold for the right ear (red mark on the hearing aid, see step 5.1) between your thumb and forefinger (*Fig. 5c*).

Take the earmold up to your ear (*Fig. 5c*). Place the ear canal part of the earmold into your ear canal (see the part circled in green (*Fig. 5c*)). Next place the hearing aid behind your ear (*Fig. 5d*). Then insert the upper part of the earmold into the upper part of the bowl of your ear (*Fig. 5e*).

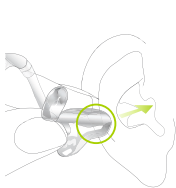


Fig. 5c



Fig. 5d



Fig. 5e

If you have problems inserting the hearing aid, use your other hand to gently pull your earlobe downwards. This opens the ear canal a little more and you can rotate the earmold until it fits correctly.

Check the fit by running your finger over the bowl of the ear. You know it fits correctly if, when doing so, you feel the contour of the bowl of the ear and not the hearing aid (*Fig. 5f*). In the beginning, please use a mirror to check.



Fig. 5f

You have now successfully inserted the right hearing aid. Now insert your left hearing aid in the same way (marked in blue, see step 5.1).

4. Using the hearing aid

Step 5.3 Inserting an UltraPower or SuperPower model with Power SlimTube with classic earmold

Follow the instructions in Step 5.2, as this procedure is identical.

Step 5.4 Inserting a CRT model with a cShell

Before inserting the cShell, check that it is correctly attached to the tube. With your right hand, place the right hearing aid (red mark on the hearing aid, see step 5.1) behind your right ear (*Fig. 5g*). Place the hearing tube between your thumb and forefinger. Insert the cShell deep enough into the ear canal so that the hearing tube lies flat on your ear (*Fig. 5h*).



Fig. 5g

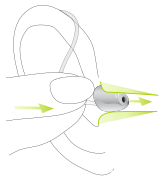


Fig. 5h

If you have problems inserting the hearing aid, use your other hand to gently pull your earlobe downwards. This opens the ear canal a little more and you can rotate the dome until it fits correctly.

You have now successfully inserted the right hearing aid. Now insert your left hearing aid in the same way (marked in blue, see step 5.1).

Step 5.5 Inserting a CRT model with a dome

Follow the instructions in Step 5.4, as this procedure is identical. If your hearing aid has a retention cord, then press it into the bowl of the ear in order to retain the dome in the ear canal (*Fig. 5i*).

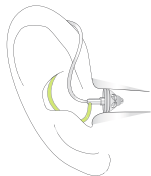


Fig. 5i

4. Using the hearing aid

Step 6. Adjusting the volume control behind the ear

Step 6.1 Adjusting the volume control behind the ear with an UltraPower or SuperPower model

To increase the volume press the volume control upwards. To decrease it, press the volume control downwards using your forefinger. Use your thumb to lightly support the hearing aid (*Fig. 6a*).

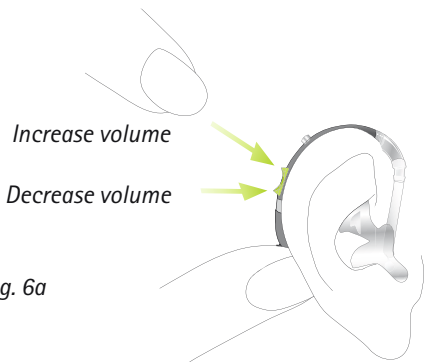


Fig. 6a

Step 6.2 Adjusting the volume control behind the ear with a CRT model

Your hearing care professional can configure your hearing aid push button as either a volume control or a program selector.

To increase the volume, press the button on your right hearing aid. To decrease the volume, press the button on your left hearing aid. Use your thumb to lightly support the hearing aid (*Fig. 6b*).

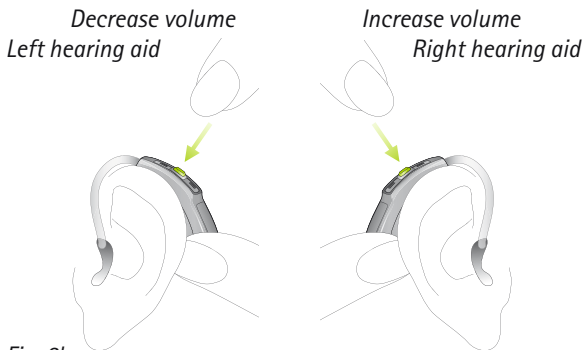


Fig. 6b

4. Using the hearing aid

Step 7. Selecting the hearing program behind the ear

Step 7.1 Selecting the hearing program behind the ear with an UltraPower or SuperPower model

In order to select or change a hearing program, press the program button displayed (*Fig. 7a*).



Fig. 7a

If you press it again, you change to the next program. A signal indicates the change.

A signal tone can only be heard if you have turned on the hearing aid and you are wearing it correctly behind the ear.

Step 7.2 Selecting the hearing program behind the ear with a CRT model

Your hearing care professional can configure your hearing aid (push) button as either a volume control or a program selector.

In order to select or change a hearing program, press the program button displayed (*Fig. 7b*).



Fig. 7b

If you press it again, you change to the next program. A signal indicates the change.

A signal tone can only be easily heard if you have turned on the hearing aid and you are wearing it correctly behind the ear.

4. Using the hearing aid

Step 8. Removing the hearing aid

Step 8.1 Removing an UltraPower or SuperPower model with classic earmold

Take the hearing aid behind your ear (do not touch the tube!) and lift it above the upper part of your ear and then gently let go (*Fig. 8a*).

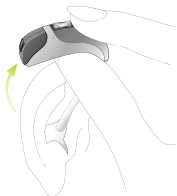


Fig. 8a

Use your thumb and forefinger to grasp the earmold in your ear (not the tube!) and gently remove it from your ear (*Fig. 8b*).

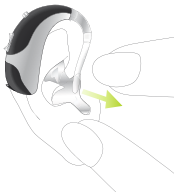


Fig. 8b

Step 8.2 Removing an UltraPower or SuperPower model with Power SlimTube with classic earmold

Follow the instructions in Step 8.1, as this procedure is identical.

Step 8.3 Removing a CRT model with a cShell

Grasp the hearing aid by the hearing tube and carefully remove it from your ear (*Fig. 8c*).

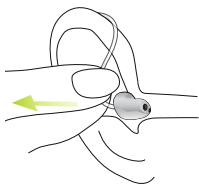


Fig. 8c

4. Using the hearing aid

Step 8.4 Removing a CRT model with dome

Follow the instructions in Step 8.3, as this procedure is identical.

⚠ In very rare cases, the dome can remain in your ear canal when removing the hearing tube from the ear. In the unlikely case that the dome does get stuck in your ear canal, it is strongly recommended to see a physician for safe removal.


Step 9. Turning off the hearing aid

You can turn off your hearing aid by opening the battery compartment (*Fig. 9a/9b*).

UP and SP model



CRT model

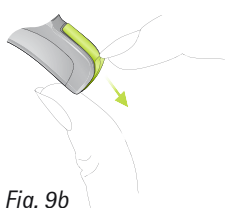

Step 10. Changing the battery

Your hearing aid will indicate with a double beep, that the battery soon needs to be changed. You typically have up to 30 minutes to change the battery, however this can differ and might be considerably shorter. We recommend that you always keep a spare battery on hand.

- ① Contrary to customary practice, particularly for pediatric consumers and for wearers of UltraPower and SuperPower products, the use of battery testers is not recommended. The battery tester may indicate that the battery is still viable while the voltage level has dropped below what the hearing aid requires for consistent function.
- ① For optimal use of the batteries, we recommend you to test and estimate the reserve time individually when you start using your hearing aid.

4. Using the hearing aid

Step 11. Changing the battery in hearing aids with a secured battery compartment

The UltraPower and SuperPower hearing aid models can be fitted with a secured battery compartment if needed for safety reasons.

Open the battery compartment by inserting a pen tip into the notch shown (*Fig. 11a*). Push the latch upwards (*Fig. 11b*). The battery compartment can now be opened.

ⓘ Do not lift the latch above 90°.

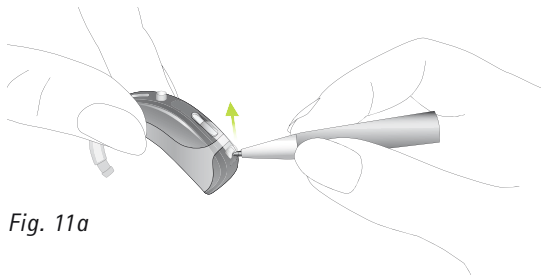


Fig. 11a

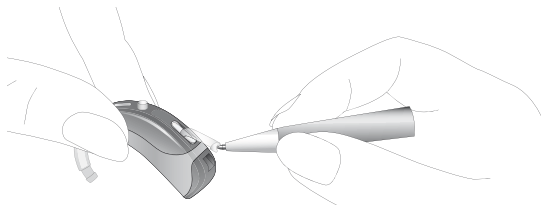


Fig. 11b

4. Using the hearing aid

Use your thumb and forefinger and pull the battery compartment to open it up (*Fig. 11c*).

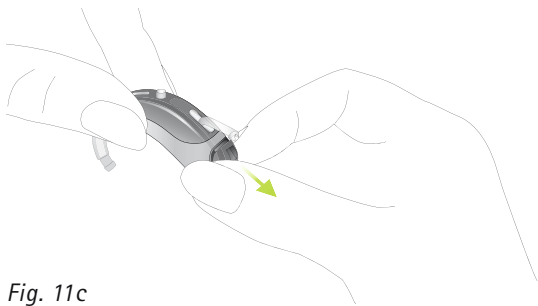


Fig. 11c

Pull the white tab with your fingernail (*Fig. 11d*) and push it in an upward direction (*Fig. 11e*). The battery can only be inserted when the tab is in upwards position (*Fig. 11e*).

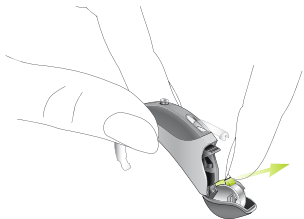


Fig. 11d



Fig. 11e

The old battery can now be replaced by a new one (*Fig. 11f*).



Fig. 11f

4. Using the hearing aid

Push the white tab down over the battery (*Fig. 11g*) so that it lies tight over the battery. Lock the battery by pushing the white tab towards the hearing aid with your finger (*Fig. 11h*). Before you close the battery compartment, check that the white tab is locked in position and the battery is secured.



Fig. 11g



Fig. 11h

"click"

Switch on the hearing aid by closing the battery compartment (*Fig. 11i*).

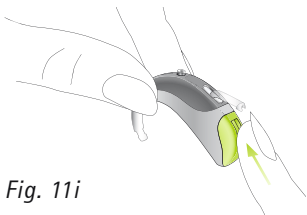


Fig. 11i

Secure the battery compartment by pushing the latch over the battery compartment (*Fig. 11j*) until it snaps into place (*Fig. 11k*). Check that the battery compartment is now secured.

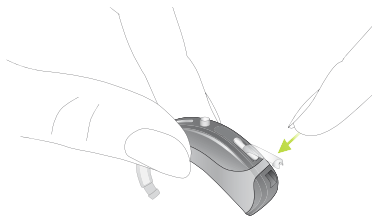


Fig. 11j

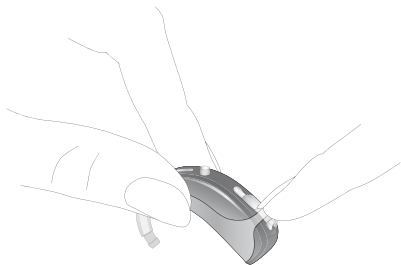


Fig. 11k

5. Care and maintenance

Diligent and routine care of your hearing aid contributes to outstanding performance and a long service life.

Please use the following specifications as a guideline:

General information

Before using hair spray or applying cosmetics, you should remove your hearing aid from your ear, because these products may damage it.

Never wash or clean the Wind and Weather Protector (see chapter 3). Doing so could cause it to lose its special acoustic features.

Naída S is WaterResistant!

Naída S hearing aids are resistant to water, sweat and dust under the following conditions:

- The battery door is fully closed. Ensure that no foreign object such as hair is caught in the battery door when it is closed.
- After exposure to water, sweat or dust, the hearing aid is cleaned and dried.

-
- The hearing aid is used and maintained as described in this user guide.
 - ① Naída S hearing aids modified with either the design-integrated FM receiver or tamperproof battery door are also WaterResistant.
 - ① Naída S hearing aids modified with a battery compartment allowing use of an audio shoe and universal FM receiver, are no longer WaterResistant.
 - ① Use of your hearing aid around water can restrict air flow to the batteries causing it to stop working. Should your hearing aid stop working after coming into contact with water, refer to the troubleshooting steps in chapter 6.

5. Care and maintenance

For an UltraPower or SuperPower model

- ① Never immerse your hearing aid in water or other liquids.
- ① Always remove your hearing aid before showering, bathing or swimming. It is not designed for any kind of swimming or diving activities.

For a CRT model

- ① Your hearing aid is resistant to temporary immersion in water, however always remove your hearing aid before diving, snorkeling or waterskiing. It is not designed for activities where it is immersed for an extended period of time or subject to excessive force.

Daily

Inspect the earmold (dome or cShell) and tube for earwax and moisture deposits. Clean the surfaces with a lint-free cloth. Never use cleaning agents such as household detergents, soap, etc. for cleaning your hearing aid. It is not recommended to rinse these parts with water, risk for lodging in the tube. If you need to clean your hearing aid intensively, ask your hearing care professional for advice and information on filters or drying capsules.

Weekly

Clean the earmold (dome or cShell) with a soft, damp cloth or with a special cleaning cloth for hearing aids. For more in depth maintenance instructions or for more than basic cleaning, please see your hearing care professional.

Monthly

Inspect your hearing tube for color changes, hardening, or cracks. In the case of such changes, the hearing tube has to be replaced. Please see your hearing care professional.

6. Troubleshooting

Hearing aids sound softer than usual/no sound

First check whether you can adjust and make the sound louder with the volume control. Then check the hearing tube and earmold for earwax residue and clean if necessary. Check whether the battery is correctly inserted and properly fitted, adjust the fit if necessary. If there is no improvement, change the battery. If the problem persists, contact your hearing care professional.

Hearing aid stops working after contact with water

■ **For an UltraPower or SuperPower model**

Open the battery compartment, dry the casing and the inside of the battery compartment with a soft cloth or tissue and replace the battery. If the problem persists, contact your hearing care professional.

■ **For a CRT model**

Rinse the hearing aid in fresh water and dry it with a soft cloth or tissue. Open the battery compartment, allow the hearing aid to dry and replace the battery. If the hearing aid uses a cShell, clean or replace the cShell's wax protection system. If the problem persists, contact your hearing care professional.

Battery discharges too quickly

Insert a new battery and make note of exactly how long it lasts. Contact your hearing care professional with this information, they will be able to help and advise you.

Crackling or buzzing noises

Check the hearing tube and earmold for earwax residue. Check the hearing tube for damage (color changes, hardening, or cracks) and contact your hearing care professional if necessary for replacement.

Whistling

Check whether the earmold is correctly and firmly fitted. If this is the case, but the problem persists, please consult with your hearing care professional.


Pain or inflammation in your ear

Remove the hearing aid and consult with your hearing care professional. If the problem is serious, please contact your physician.

7. Wireless Accessories and FM systems

7.1 Wireless accessories

Phonak has developed a broad range of wireless accessories for a convenient control of your hearing aid. It includes remote controls, wireless audio streaming devices for connection to TV (via TVLink), cell phones etc. refer to schematic below:

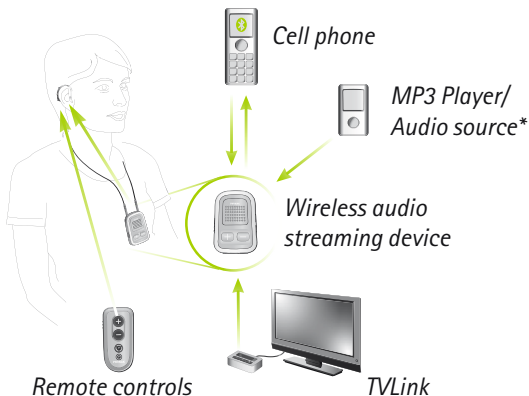

* Non wireless audio sources can be connected by wiring directly to the wireless audio streaming device.


7.2 FM systems

FM systems are wireless communication systems, consisting of a transmitter and a receiver. For more information, please visit www.phonak.com or ask your hearing care professional. If you already have an FM system and need additional information, please consult the relevant user instructions.

7. Wireless accessories and FM systems

Transmitter: placed close to the speaker or connected to the sound source (Cell phone, TV etc.), transmitting the sound over distance to the receiver directly without disruptive background noises.

Receiver: worn by the listener (neckloop) or connected directly to the Behind-The-Ear hearing aids.

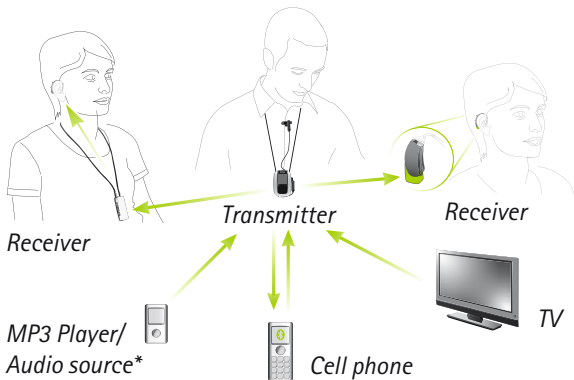

* Non wireless audio sources can be connected by wiring directly to the transmitter.


7.2.1 FM system description

The following pages describe the two FM receiver models: the **design-integrated FM receiver** and the **audio shoe with universal FM receiver**.


On the following pages use the labeled drawings to identify the FM receiver that you use.

7.2.2 Design-integrated FM receiver

UP and SP model

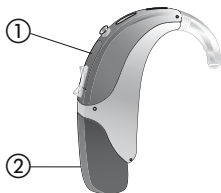

CRT model

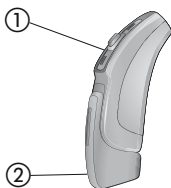

- ① Hearing aid
- ② Design-integrated FM receiver

7. Wireless accessories and FM systems


With the design-integrated FM receiver attached the hearing aid offers all functionalities and usability as without the FM receiver (see chapter 4).

- ① Your hearing care professional has already prepared your hearing aid for use with the design-integrated FM receiver. The design-integrated FM receiver can not be removed. Please ask you hearing care professional for further information.


Naída S with the design-integrated FM receiver is WaterResistant. For more information on the water resistant properties of the Naída S housing, please see chapter 5.

7.2.3 Audio shoe with universal FM receiver

UP and SP model

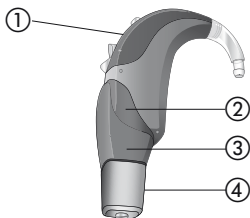

CRT model

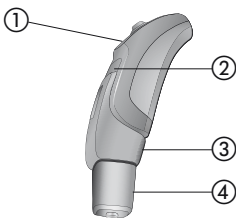

- ① Hearing aid
- ② Battery compartment
- ③ Audio shoe for universal FM receiver
- ④ Universal FM receiver

With the universal FM receiver attached the hearing aid offers all functionalities and usability as without the FM receiver (see chapter 4).

7. Wireless accessories and FM systems

- ① Naída S with the audio shoe and the universal FM receiver is not WaterResistant. Please ask you hearing care professional for further information.

Attaching the audio shoe to the hearing aid

Your hearing care professional has already prepared your hearing aid for use with a universal FM receiver. Figures 7.1 to 7.5 show the **UltraPower** and **SuperPower model**, follow these instructions for the **CRT model**, as this procedure is identical.

To attach the audio shoe, hold the hearing aid between the thumb and forefinger of one hand and the audio shoe in the other hand. Push the audio shoe towards the bottom of the hearing aid. Match the grooves of the hearing aid with those of the audio shoe (*Fig. 7.1*).

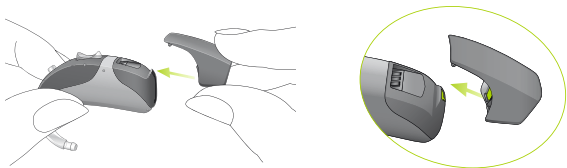


Fig. 7.1

Push the hearing aid down until it snaps into place (Fig. 7.2). The audio shoe is now properly attached to the hearing aid. A universal FM receiver can now be attached.

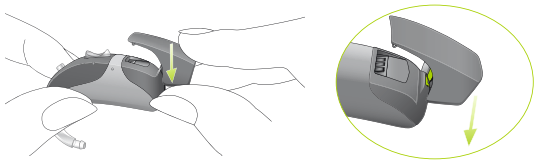


Fig. 7.2

Attaching the universal FM receiver to the audio shoe
Connect the universal FM receiver to the audio shoe by pushing the three pins of the universal FM receiver into the three sockets of the audio shoe. The universal FM receiver is correctly connected to the audio shoe once it is clicked in place (Fig. 7.3). Please refer to the specific FM receiver user instruction.

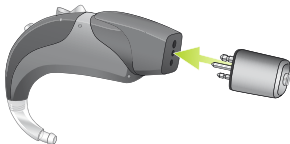


Fig. 7.3

7. Wireless accessories and FM systems

- ① To have FM immediately active when required, Phonak recommends to leave the universal FM receiver attached when not in use.

Remove the audio shoe from the hearing aid

To remove the audio shoe, take the hearing aid between your thumb and forefinger. Pull the locking system on the audio shoe back (*Fig. 7.4*).

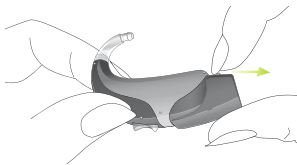


Fig. 7.4

At the same time pull the hearing aid upwards with the other hand (*Fig. 7.5*). The audio shoe is now removed from the hearing aid.

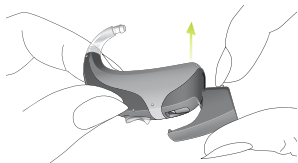


Fig. 7.5

8. Service and warranty

8.1 Local warranty

Please ask the hearing care professional, where you purchased your hearing aid, about the terms of the local warranty.

8.2 International warranty

Phonak offers a one year limited international warranty, valid starting from the date of purchase. This limited warranty covers manufacturing and material defects in the hearing aid itself, but not accessories such as batteries, tubes, earmolds, external receivers. The warranty only comes into force if a proof of purchase is shown.

The international warranty does not affect any legal rights that you might have under applicable national legislation governing sale of consumer goods.

8. Service and warranty

8.3 Warranty limitation

This warranty does not cover damage from improper handling or care, exposure to chemicals or undue stress. Damage caused by third parties or non-authorized service centers renders the warranty null and void. This warranty does not include any services performed by a hearing care professional in their office.

Serial number (left side):

Serial number (right side):

Date of purchase:

Authorized hearing care
professional (stamp/
signature):

9. Compliance information

Europe:

Declaration of Conformity

Hereby Phonak AG declares that this Phonak product meets the requirements of the Medical Devices Directive 93/42/EEC as well as the Radio and Telecommunications Terminal Equipment Directive 1999/5/EC. The full text of the Declaration of Conformity can be obtained from the manufacturer or the local Phonak representative whose address can be taken from the list on www.phonak.com (Phonak worldwide locations).

Australia/New Zealand:



R-NZ

Indicates a device's compliance with applicable Radio Spectrum Management's (RSM) and Australian Communications and Media Authority (ACMA) regulatory arrangements for the legal sale in New Zealand and Australia. The compliance label R-NZ is for radio products supplied in the New Zealand market under conformity level A1.

9. Compliance information

The wireless models listed in chapter 1 are certified under:

UP and SP models:

USA FCC ID: KWC-WHSBTE1U

Canada IC: 2262A-WHSBTE1U

CRT models:

USA FCC ID: KWC-WHSSAN1

Canada IC: 2262A-WHSSAN1

Notice 1:

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- 1) this device may not cause harmful interference, and
- 2) this device must accept any interference received, including interference that may cause undesired operation.

Notice 2:

Changes or modifications made to this device not expressly approved by Phonak may void the FCC authorization to operate this device.

Notice 3:

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This device generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the device and receiver.
- Connect the device into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

10. Information and explanation of symbols



With the CE symbol, Phonak AG confirms that this Phonak product – including accessories – meets the requirements of the Medical Devices Directive 93/42/EEC as well as the R&TTE Directive 1999/5/EC on radio and telecommunications equipment. The numbers after the CE symbol correspond to the code of certified institutions that were consulted under the above-mentioned directives.



This symbol indicates that the products described in these user instructions adhere to the requirements for an applied part of Type B of EN 60601-1. The surface of the hearing aid is specified as an applied part of Type B.



During transportation keep dry.



This symbol indicates that it is important for the user to read and take into account the relevant information in these user guides.



This symbol indicates that it is important for the user to pay attention to the relevant warning notices in these user guides.



Important information for handling and product safety.

Operating conditions

The product is designed such that it functions without problems or restrictions if used as intended, unless otherwise noted in these user guides.

10. Information and explanation of symbols

Transport and storage conditions

Temperature: -20° to $+60^{\circ}$ Celsius
(-4° to $+140^{\circ}$ Fahrenheit)

Humidity transportation: Up to 90%
(non condensing)

Humidity storage: 0% to 70%, if not
in use. See instruction in chapter 2.2
regarding drying the hearing aid after
use.

Atmospheric pressure: 200 hPa to
1500 hPa



The symbol with the crossed-out garbage bin is to make you aware that this hearing aid may not be thrown away as normal household waste. Please dispose of old or unused hearing aids, at waste disposal sites intended for electronic waste, or give your hearing aid to your hearing care professional for disposal. Proper disposal protects the environment and health.



Notes



Manufacturer:

Phonak AG

Laubisrütistrasse 28

CH-8712 Stäfa

Switzerland

www.phonak.com

CE
0459



029-1061-02/M6.00/2016-02/A+W/FO Printed in Switzerland
© Phonak AG All rights reserved